

## HOW TO USE

The medication needs multiple layers on the lesion site: it can, therefore, be folded several times. Subsequently it requires a secondary dressing (such as, for instance, non-woven fabric, silver-free calcium alginates or hydrofibres, if particularly exuding) which does not absorb its oily component.

1. Cleanse the lesion site with saline and dab dry with sterile gauze.
2. Apply the medication on the lesion and adapt as needed.

Please, read the instructions and warnings for use carefully before using the device.


## AVAILABLE SIZES

### Cotton and polyamide elastic bandage

- NBE4.5x8S** Pack of 1 piece 4,5 m x 8 cm
- NBE4.5x8** Pack of 10 pieces 4,5 m x 8 cm

### Bibliography

- Montecchi et al., New Microbiologica 2013; 36: 289-302.
- Cassino R., et al., MINERVA CHIR 2015; 70: 23-31
- Laureti S., et al., Tech Coloproctol. 2016; 20: 139-40.
- Vidotto G., et al. Gastrointestinal Nursing 2020; Sup4: S23-S30
- Esposito C., et. al., JLAST 2020; 4: 1-7
- Esposito C., et. al., Front. Pediatrics 2020; 8: 345(1-9)

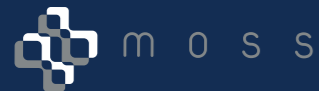

**MOSS S.p.A.** Via all'Erno, 5 - 28040 Lesa (NO) - Italy  
Tel. +39 0322 772032 - Fax +39 0322 292733  
moss@moss-brand.it - www.moss-brand.com

### Distributor in Ireland and Northern Ireland:



Unit 31 Southern Cross Business Park  
Boghall Road, Bray, Co. Wicklow, A98 PX00  
ROI +353 1 201 4620 - NI +44 28 9099 5160  
sales@fleetwoodhealthcare.ie - www.fleetwoodhealthcare.ie





Class IIb medical device.

Dressing in the form of bandage for the treatment of skin lesions, impregnated with an oleic matrix-based gel releasing **Reactive Oxygen Species (ROS)**.

## MECHANISM OF ACTION

- 1 **Novox®-Roll** is a medical device used in the treatment of ulcerative or traumatic skin lesions, that promotes the wound healing process by creating a micro-environment favourable to the activation of microcirculation.
- 2 **Novox®-Roll** is able to build a protective barrier and a micro-environment unfavourable to the possible proliferation of pathogens.
- 3 **Novox®-Roll** has a film-forming and soothing effect due to the oily nature of the product, combined with the healing and re-epithelizing properties of oxygen-enriched oil.
- 4 Due to its specific conformation, **Novox®-Roll** allows to treat large-size lesions.

## TECHNOLOGY

**Novox®-Roll** is based on a technological process studied and perfected by the R&D team of MOSS, which has maximized the enrichment of oxygen as part of an oleic matrix, so as to lend the medical device its technical features of constant, continuous release of **ROS**.

## INDICATIONS

- 1 Skin sampling areas
- 2 Leg ulcers
- 3 Pressure ulcers
- 4 Traumatic or stab wounds
- 5 Post-surgical wounds
- 6 Fistulas
- 7 Burns

## ADVANTAGES

- 1 Free from active pharmaceutical ingredients (APIs) and preservatives
- 2 Does not induce antibiotic resistance
- 3 Uniform distribution over the wound area, facilitated by the flexibility of the product
- 4 Simplified application on large-size sleeve lesions or wounds involving interdigital areas of the foot and part of the limb
- 5 Promotes the drainage of exudate and protects surface lesions from adhesions
- 6 Can be packed in cavity lesions and removed with forceps
- 7 Can be cut out and shaped according to the conditions of the lesions
- 8 If desired, can be used in association with compression bandages
- 9 Affords atraumatic removal