



## HOW TO USE

1. Cleanse the lesion with saline according to the protocol in use and dab dry.
2. Apply the medication on the lesion, to obtain a uniform layer a few millimeters thick.
3. Depending on the condition of the lesion, a greasy or non-adherent gauze is recommended as secondary dressing, to hold the medication in place on the wound bed.

Please, read the instructions and warnings for use carefully before using the device.

## AVAILABLE SIZES

- NS005** Pack of 10 x 5 ml syringes  
**NS005S** Pack of 1 x 5 ml syringe

## Bibliography

- Montevecchi et al., New Microbiologica 2013; 36: 289-302.
- Cassino R., et al., MINERVA CHIR 2015; 70: 23-31
- Laureti S., et al., Tech Coloproctol. 2016; 20: 139-40.
- Vidotto G., et al. Gastrointestinal Nursing 2020; Sup4: S23-S30
- Esposito C., et. al., JLAST 2020; 4: 1-7
- Esposito C., et. al., Front. Pediatrics 2020; 8: 345(1-9)

n o v o x<sup>®</sup>



**MOSS S.p.A.** Via all'Erno, 5 - 28040 Lesa (NO) - Italy  
Tel. +39 0322 772032 - Fax +39 0322 292733  
moss@moss-brand.it - www.moss-brand.com

## Distributor in Ireland and Northern Ireland:




Unit 31 Southern Cross Business Park  
Boghall Road, Bray, Co. Wicklow, A98 PX00  
ROI +353 1 201 4620 - NI +44 28 9099 5160  
sales@fleetwoodhealthcare.ie - www.fleetwoodhealthcare.ie

CE 0477  
Rev. 10-IR May/2021

n o v o x<sup>®</sup>

Class IIb dressing for the treatment of skin lesions through the extended release of **Reactive Oxygen Species (ROS)**.

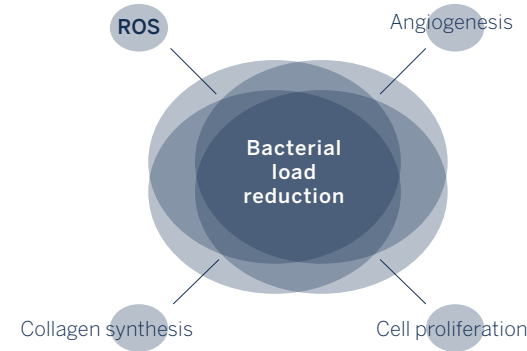

## MECHANISM OF ACTION

- 1 **Novox<sup>®</sup>** is a medical device used in the treatment of ulcerative or traumatic skin lesions, that promotes the wound healing process by creating a micro-environment favourable to the activation of microcirculation.
- 2 **Novox<sup>®</sup>** is able to build a protective barrier and a micro-environment unfavourable to the possible proliferation of pathogens.
- 3 **Novox<sup>®</sup>** has a film-forming and soothing effect due to the oily nature of the product, combined with the healing and re-epithelizing properties of oxygen-enriched oil.

## TECHNOLOGY

**Novox<sup>®</sup>** is based on a technological process studied and perfected by the R&D team of MOSS, which has maximized the enrichment of oxygen as part of an oleic matrix, so as to lend the medical device its technical features of constant, continuous release of **ROS**.

## INDICATIONS

- 1 Traumatic ulcers
- 2 Vascular ulcers
- 3 Pressure ulcers
- 4 Fistulas and post-surgical non-healing wounds
- 5 Traumatic wounds (such as, for instance, abrasions and lacerations)
- 6 Skin and mucosal lesions

## ADVANTAGES

- 1 Free from active pharmaceutical ingredients (APIs) and preservatives
- 2 Does not induce antibiotic resistance
- 3 Simplified application and removal
- 4 Uniform distribution over the wound bed, facilitated by the consistency of the product
- 5 Targeted application by means of the dedicated spout

## EFFECTIVENESS AGAINST MICROBIAL CONTAMINATION

